



Mouse TGF- β 1 ELISA

Catalog Number EA-2521

(For Research Use Only)

Introduction

The transforming growth factor beta 1 (TGF- β 1) gene codes a multifunctional cytokine that controls proliferation, differentiation, and other functions in many cell types, including cancer cells, the surrounding stromal cells, immune cells, endothelial and smooth-muscle cells. It causes immunosuppression and angiogenesis, which makes the cancer more invasive. TGF- β also converts effector T-cells, which normally attack cancer with an inflammatory (immune) reaction, into regulatory (suppressor) T-cells, which turn off the inflammatory reaction. TGF- β induces apoptosis in numerous cell types. TGF- β can act on adipocyte precursor cells (1). TGF- β 1 has been shown to be a potent inhibitor of the differentiation of adipogenic cell lines (2). In addition, a differentiation-defective, insulin-independent cell line derived from the adipogenic cell line 1246 produces in its conditional medium a TGF- β 1-like polypeptide which could modulate the cell ability to differentiate in an autocrine fashion. Increased TGF- β 1 expression was associated with BMI and abdominal adipose tissue in morbid obesity (4).

Principle of the assay

TGF- β 1 ELISA is based on the principle of a solid phase enzyme-linked immunosorbent assay. The assay utilizes a mouse anti-mouse TGF- β 1 antibody for immobilization on the microtiter wells and chicken anti-mouse TGF- β 1 antibodies along with streptavidin conjugated to horseradish peroxidase (HRP) for detection. The test sample is allowed to react simultaneously with the two antibodies, resulting in the TGF- β 1 molecules being sandwiched between the solid phase and enzyme-linked antibodies. After incubation, the wells are washed to remove unbound-labeled antibodies. A HRP substrate, TMB, is added to result in the development of a blue color. The color development is then stopped with the addition of Stop Solution changing the color to yellow. The concentration of TGF- β 1 is directly proportional to the color intensity of the test sample. Absorbance is measured spectrophotometrically at 450 nm.



Diagram of ELISA

Materials provided with the kit

- 96 well microplate coated with a mouse anti-mouse TGF- β 1 antibody (4°C).
- Biotin labeled chicken anti-mouse TGF- β 1 antibodies (-20°C).
- Streptavidin-HRP conjugate (4°C).
- Recombinant TGF- β 1 standard (-20°C).
- 1X Diluent buffer (4°C).
- 5X Assay wash buffer (RT)
- Substrate (4°C).
- Stop Solution (4°C).

Material required but not provided

- Microplate reader capable of measuring absorbance at 450 nm
- Deionized or distilled water.


Reagent preparation before starting experiment

- Dilute the 5x Assay wash buffer to 1x buffer
40ml 5x Assay wash buffer
160ml ddH₂O
- Dilute 50 times of mouse recombinant TGF- β 1 (220ng/ml) with 1X Diluent buffer to 4400pg/ml and then 2-fold serial dilutions. Dilute 50 times by adding 4ul Mouse Recombinant TGF- β 1 in 200ul 1X Diluent Buffer (See Step 2 in "Assay Procedure" for detailed procedure)
- Dilute 400 times of biotin labeled chicken anti-mouse TGF- β 1 antibody with 1X Diluent buffer before use.
- Dilute 200 times of streptavidin-HRP with 1X Diluent buffer before use.

7. Add 100 μ l of diluted streptavidin-HRP conjugate to each well and incubate for 45 min at room temperature with gentle shaking.
8. Repeat the aspiration/wash as in step 4.
9. Add 100 μ l of substrate to each well and incubate for 10-30 minutes.
10. Add 50 μ l of Stop solution to each well. The color in the wells should change from blue to yellow.
11. Determine the optical density of each well with a microplate reader at 450 nm within 30 minutes.

Assay procedure

1. Cut the sealing film over the plate and remove it from the desired number of well strips. Make sure the rest of wells are well sealed.
2. See instruction and diagram below for standard preparation.

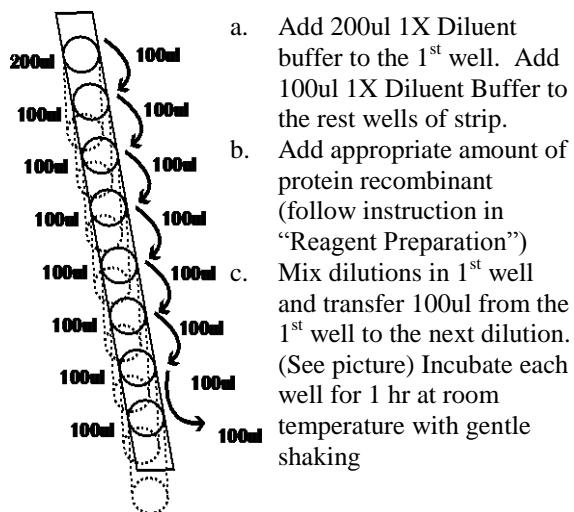

3. Add 100ul of sample per well and incubate for 1 hour at room temperature with gentle shaking.
4. Aspirate each well and wash by adding 200 μ l of 1X Assay wash buffer. Repeat the process three times for a total of three washes. Complete removal of liquid at each wash. After the last wash, remove any remaining liquid by inverting the plate against clean paper towels.
5. Add 100 μ l of diluted biotin-labeled anti-mouse TGF- β 1 antibody to each well and incubate for 1 hour at room temperature with gentle shaking.
6. Repeat the aspiration/wash as in step 4.