



Human Anti-SSA (Ro-52) ELISA Kit

Catalog Number EA-5004

(For Research Use Only)

Introduction

Antibodies to the SSA antigen (also known as Ro antigen) are one of the most frequent serological markers of autoimmunity in rheumatic diseases. They appear in 60-70% of patients with Sjögren's syndrome (SS), 30-40% of patients with systemic lupus erythematosus (SLE), and 3-5% of patients with rheumatoid arthritis (RA). The SSA antigen is comprised of an acidic 52 kDa protein that may exist in complex with RNA (80-112 bases). However, RNA is not required for SSA antigenicity. The SSA antigen is predominantly cytoplasmic, and not located in the nucleus. Therefore, patients with antibody to SSA may be ANA negative on routine testing. Approximately 62% of ANA negative lupus patients have antibodies to SSA. Two types of anti-Ro/SSA antibodies have been identified; Ro/SSA antigens of 60 kDa and 52 kDa. Anti-SSA-52 kDa antibodies are linked to certain disorders such as SS, SLE, neonatal lupus and congenital heart block. Clinically, the presence of SSA52 has been reported in a wide variety of diseases, includes inflammatory myositis, primary biliary cirrhosis and SS.

Principle of the assay

Anti-SSA ELISA kit measures anti-SSA antibodies in the serum. It is based on the principle of a solid phase enzyme-linked immunosorbent assay. The assay utilizes SSA protein for immobilization on the microtiter wells and anti-human IgG antibodies conjugated to horseradish peroxidase (HRP) for detection. The test sample is allowed to react simultaneously with the two components, resulting in anti-SSA antibodies being sandwiched between the solid phase and enzyme-linked antibodies. After incubation, the wells are washed to remove unbound-labeled antibodies. A HRP substrate, TMB, is added to result in the development of a blue color. The color development is then stopped with the addition of Stop Solution changing the color to yellow. The concentration of anti-SSA is directly proportional to the color intensity of the test sample. Absorbance is measured spectrophotometrically at 450 nm.

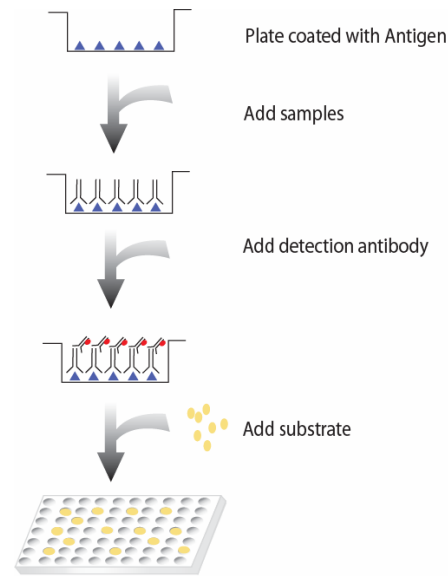


Diagram of ELISA

Materials provided with the kit

- 96-well plate coated with SSA (Ro-52) (4°C).
- Anti-human IgG antibody conjugated to HRP (4°C).
- 1X Diluent buffer (4°C).
- 5X Assay wash buffer (4°C).
- Substrate (4°C).
- Stop Solution (4°C)

Material required but not provided

- Microplate reader capable of measuring absorbance at 450 nm
- Shaker

Reagent preparation before starting experiment

- Dilute the 5x Assay wash buffer to 1x buffer
40ml 5x Assay wash buffer
160ml ddH₂O
- Dilute 1000 times of anti-human IgG antibody conjugated to HRP with 1X Diluent buffer.

Storage and Preparation

Store all reagents at 2-8 °C.

All reagents must be brought to room temperature (20-25 °C) prior to use.

When stored at 2-8 °C, the diluted Assay wash buffer is stable until the kit expiration date.

Precautions

Human blood derivatives and patient specimens should be considered potentially infectious. All human derived components need to be tested for the negative HBsAg, HCV, HIV-1 and 2 and HTLV-I. Follow good laboratory practices in storing, dispensing and disposing of these materials.

Assay procedure

1. Cut the sealing film over the plate and remove it from the desired number of well strips. Make sure the rest of wells are well sealed.
2. Add 100 µl of diluted samples or positive control (1:100 diluted with 1X Diluent Buffer) per well and incubate for 1 hour at room temperature with gentle shaking.
3. Aspirate each well and wash by adding 200µl of 1X Assay wash buffer. Repeat the process twice for a total of three washes. Completely remove liquid at each wash by firmly tapping the plate against clean paper towels.
4. Add 100µl of diluted anti-human IgG antibody conjugated to HRP to each well and incubate for 0.5 hours at room temperature with gentle shaking.
5. Repeat the aspiration/wash as in step 3.
6. Add 100 µl of Substrate to each well and incubate for 5-30 minutes.
7. Add 50µl of Stop solution to each well. The color in the wells should change from blue to yellow.
8. Determine the optical density of each well with a microplate reader at 450 nm within 30 minutes.

Example

